



# CORYNEBACTERIUM MINUTISSIMUM

## qPCR TEST HANDBOOK

For Research Use Only. Not intended for diagnostic use.

CAT NO.: YS-qP-EC-C.minut-100

100 reactions  
with Endogenous Control and Frozen MasterMix

VERSION 8.3



**YouSeq Ltd**  
8 Moorside Place  
Moorside Road  
Winchester  
SO23 7FX  
United Kingdom

**+44 (0) 333 577 6697**  
**hello@youseq.com**

**youseq.com**



## INTENDED USE

This qPCR test offers an efficient and user-friendly solution for the detection of *Corynebacterium minutissimum* (C.minut) in extracted nucleic acid samples from a variety of sources. It is intended for use by trained professionals in a suitable molecular biology laboratory.

## SPECIFICITY AND SENSITIVITY

### Specificity


The YouSeq qPCR test for detection of *Corynebacterium minutissimum* is designed to have the broadest detection profile possible and detect all clinically relevant strains. The primers and probes typically have a  $\geq 95\%$  homology with all reference data used, from relevant, publicly available databases at the time of design.

The target gene for *Corynebacterium minutissimum* (gyrB gene) has been demonstrated to have a distinctive sequence making it an ideal target for highly specific detection.

For further information on the detection profile of the product, please do not hesitate to contact our team: [support@youseq.com](mailto:support@youseq.com)

### Sensitivity

The qPCR test is suitable for the detection of *Corynebacterium minutissimum*, across a wide dynamic range. Under ideal PCR conditions the assay can detect less than 100 copies of the target in the PCR reaction.

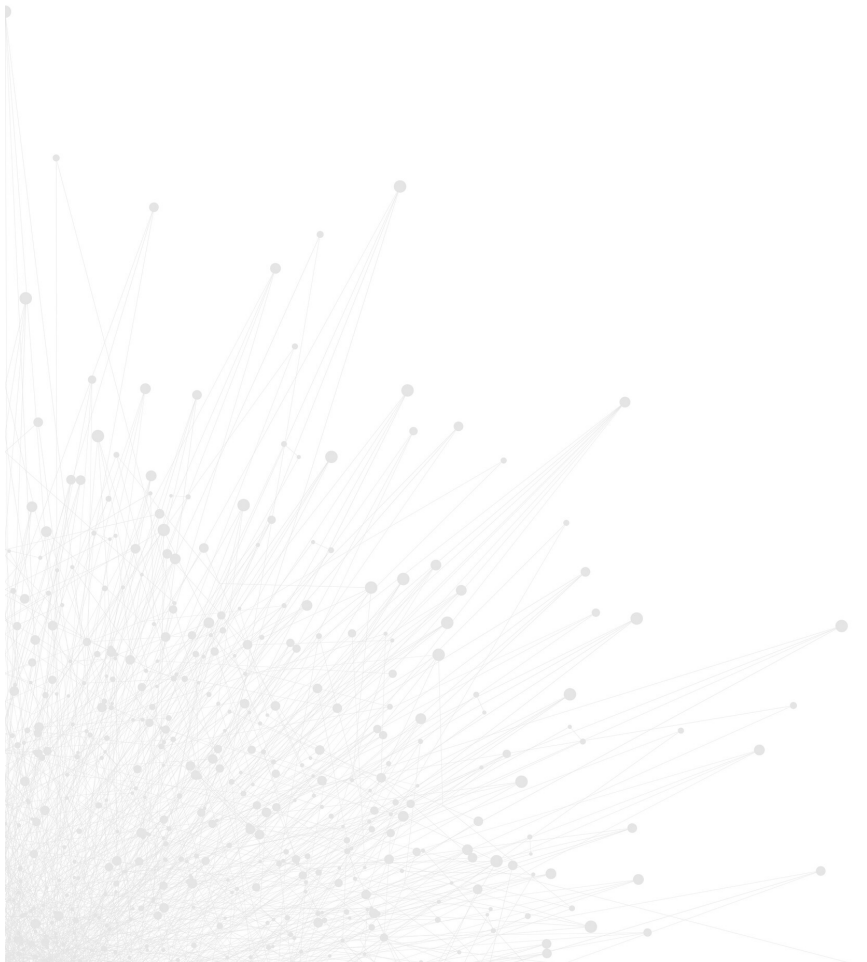



YOUSEQ

## CONTENTS

| Component                                                    | Cap Colour                                                                          | Volume  |
|--------------------------------------------------------------|-------------------------------------------------------------------------------------|---------|
| Corynebacterium minutissimum specific primer/probe mix (FAM) |  | 100 µL  |
| Endogenous Control primer/probe mix (VIC/HEX)                |  | 100 µL  |
| Tetra™ 2X qPCR MasterMix                                     |  | 1 mL    |
| Corynebacterium minutissimum Positive Control (PTC)          |  | 500 µL* |
| Template Resuspension Buffer (TRB)                           |  | 1.5 mL  |
| DNase/RNase Free Water                                       |  | 1.5 mL  |
| ROX Passive Reference Dye                                    |  | 10 µL   |

\* Supplied dried – requires resuspension. See instructions in resuspension section.

## RECOMMENDED ADDITIONAL REAGENTS & MATERIALS

Nucleic Acid extraction kit

General laboratory equipment (pipettes, pipette tips, (micro)centrifuge tubes, compatible strip tubes/plates, plate seals, etc.).

qPCR instrument with channels to detect FAM and VIC/HEX.

## BEST PRACTICE

### Decontamination:

Before beginning laboratory work, thoroughly decontaminate any work surfaces and pipettes being used, to eliminate potential contamination.

### General use and set-up:

All components should be fully defrosted with contents at the bottom of the tube before opening. To ensure contents are at the bottom, centrifuge or gently tap the tube. After use, reagents should be returned to the freezer.

Once any reagents are resuspended, mark the tick box on the tube for future reference. After this, or after combining reagents, the tube should be pulse vortexed to ensure it is mixed well.

It is advised to set up the tubes/plate and reaction mix on ice to minimise artefact formation, which may reduce sensitivity.

When preparing the qPCR reaction mix, it is recommended to incorporate an overage when calculating the total number of reactions to compensate for potential volume losses incurred during pipetting.

### Set-up environments:

It is best practice to set up qPCR tubes/plates in two different environments - a clean (no template) lab and PCR (template) lab.

No Template Control(s) (NTC) and Positive Control(s) (PTC) should be included in every run. To reduce contamination, NTCs and samples can be set up and sealed in a clean lab before moving to the PCR lab.



YOUSEQ

FI

LI

FE

LE

BDE

BRE

BDI

BRI

## BENCH SIDE PROTOCOL

### RESUSPENSION

Before first use, resuspend the designated component with the correct reagent and specified volume, as per the table below:

1. Add the resuspension reagent and pulse vortex the tube to ensure each is mixed well.

| Component                                           | Reagent | Volume | Location |
|-----------------------------------------------------|---------|--------|----------|
| Corynebacterium minutissimum Positive Control (PTC) | TRB     | 500 µL | PCR lab  |

### ROX (INSTRUMENT DEPENDENT)

ROX is required for instruments that use ROX as a passive reference. The table below outlines the qPCR instruments that require the addition of ROX.

1. If ROX is required:
  - a. Dilute the ROX supplied according to the table below based on the intended qPCR instrument to be used.
  - b. Directly transfer 5µL of prepared ROX to the Tetra™ MasterMix.

| Level of ROX      | Instruments                                                            | Step 1: Volume of water to add to ROX tube | Step 2: Add to MasterMix tube |
|-------------------|------------------------------------------------------------------------|--------------------------------------------|-------------------------------|
| High ROX:         | Applied Biosystems 7000, 7300, 7700, 7900, StepOne, StepOne Plus       | No Dilution Required                       | 5 µL                          |
| Low ROX:          | Applied Biosystems 7500 & 7500 FAST, ViiA7, Quantstudio, Stratagene MX | 130 µL                                     | 5 µL                          |
| ROX Not Required: | All Other Instruments                                                  | Not Required                               | Not Required                  |

## qPCR REACTION SET-UP

1. Retrieve the required components and appropriate plasticware for qPCR reaction set-up.
2. In an appropriately sized (micro)centrifuge tube, combine the following reagents to create a reaction mix that will cover all required tube/wells (e.g. **samples**, **NTC** and **PTC**).

**Please note:** When calculating required reactions, include an overage to allow for volume loss during pipetting

| Component                                              | Volume (per reaction) |
|--------------------------------------------------------|-----------------------|
| Tetra™ 2X qPCR MasterMix                               | 10 µL                 |
| Corynebacterium minutissimum specific primer/probe mix | 1 µL                  |
| Endogenous Control primer/probe mix                    | 1 µL                  |
| Reaction mix volume                                    | 12 µL                 |

3. Mix the combined reagents by briefly vortexing or inverting.
4. Dispense 12 µL of the reaction mix into all required tubes/wells.
5. For the **NTC(s)**, add 8 µL of DNase/RNase Free Water into required tube/well(s).
6. For each extracted **sample**, add 8 µL into required tube/well(s).

**Please note:** It is best practice to seal the **NTC(s)** and **sample** tubes/wells before proceeding to the positive control steps.


7. For the **positive controls(s)**, follow either step a. or b. below:
  - a. For **qualitative** results, add 5 µL of the resuspended **PTC** (tube 1) and 3 µL of DNase/RNase Free Water, into each designated positive control tube/well. Proceed to step 8.
  - b. For **quantitative** results, perform a serial dilution of the resuspended **PTC** (tube 1) to create a six-point standard curve:
    - I. Add 45 µL of Template Resuspension Buffer into 5 fresh microcentrifuge tubes and label them 2, 3, 4, 5 and 6.
    - II. Pipette 5 µL of **PTC** into tube 2.
    - III. Mix by pipetting up and down 10 times.
    - IV. Change pipette tip and pipette 5 µL from tube 2 into tube 3.
    - V. Mix by pipetting up and down 10 times.
    - VI. Repeat steps 'IV' and 'V' to complete the dilution process for tubes 4, 5 and 6.
    - VII. Add 3 µL of DNase/RNase Free Water into each designated positive control tube/well.
    - VIII. Add 5 µL of the required **PTC** dilution (tubes 1 – 6), into each designated positive control tube/well.

**Please note:** The described standard curve provides a dynamic range as outlined in the table below.

| Tube No. | Copies of Target/reaction |
|----------|---------------------------|
| 1        | 1,000,000                 |
| 2        | 100,000                   |
| 3        | 10,000                    |
| 4        | 1,000                     |
| 5        | 100                       |
| 6        | 10                        |

8. Seal the tube(s)/plate.

9. Briefly vortex the sealed tube(s)/plate, then spin in a centrifuge to ensure all reagents are fully resuspended and at the bottom of the tubes/wells before proceeding.

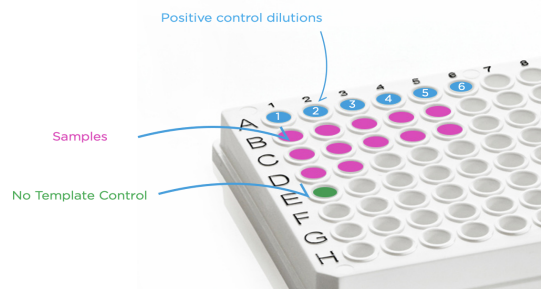

## qPCR AMPLIFICATION PROTOCOL

1. Load the plate onto the qPCR instrument and set up the qPCR protocol following the table below.
2. Set the total reaction volume to 20  $\mu$ L.

| Temperature                                                                            | Time       | Number of Cycles |
|----------------------------------------------------------------------------------------|------------|------------------|
| 95°C                                                                                   | 3 minutes  | -                |
| 95°C                                                                                   | 15 seconds | x 45             |
| 60°C  | 60 seconds |                  |

 Collect fluorogenic data through FAM and VIC/HEX channels during this step.

3. Start the run.

## INTERPRETATION OF RESULTS - OVERVIEW

If using single threshold analysis – YouSeq recommends setting the threshold at 10% of the End Point Fluorescence (EPF) for each channel.

- For the Target channel, use the **PTC** (tube 1) EPF to set the threshold.
- For the Endogenous Control, use the average EPF from **samples** to set the threshold.

### Results interpretation:

| Reaction Type             |                               | qPCR Signal  |              |
|---------------------------|-------------------------------|--------------|--------------|
| Positive control (tube 1) |                               | 18.5 $\pm$ 2 | 18.5 $\pm$ 2 |
| No template control       |                               | -            | -            |
| Sample                    | Endogenous Control: (VIC/HEX) | $\leq$ 28    | $\leq$ 28    |
|                           | Target: (FAM)                 | +            | -            |

Result

Positive result

Negative result

## INTERPRETATION OF RESULTS – CONTROLS EXPLAINED

### Positive control

The **PTC** (tube 1) should amplify in a Cq range of  $18.5 \pm 2$  for the target. If this Cq range is not achieved, the sample test result is invalid and should be repeated.

To achieve the most accurate quantitative result from the positive control standard curve, an efficiency between 90% to 110% is desirable. If it falls outside of this range, preparing a fresh standard curve and repeating the run may improve the efficiency value.

**Please note:** The positive control is a sequence representative of the target region and does not contain the organism's entire genome. The positive control does not include the Endogenous Control sequence and should not be expected to amplify in the VIC/HEX channel.

### No template control

The **NTC** tube/well(s) should be negative, with no amplification.


**Please note:** Background laboratory contamination can result in a very late signal in **NTC** tube/well(s). If the **NTC** has amplification, comparison to the sample test result is necessary:

- If the **NTC** is  $\geq 5$  Cq later than the sample signal, it can be considered negative, and the sample test result is valid.
- If the **NTC** is  $< 5$  Cq later, the sample test result is invalid – the test should be repeated after potential sources of contamination have been removed.

### Endogenous Control

Detection of the Endogenous Control is through the VIC/HEX channel. The Cq value from the Endogenous Control will vary according to the amount of DNA in the sample. Cq values  $\leq 28$  indicate a successful extraction has taken place with sufficient host derived DNA in the sample. If the signal is later than this, repeating the nucleic acid extraction is advised.

**Please note:** If the sample delivers a strong positive result for the target of interest, then amplification of the Endogenous Control may be affected and may appear later. In this case, lack of Endogenous Control amplification is acceptable, and the sample test result is valid.

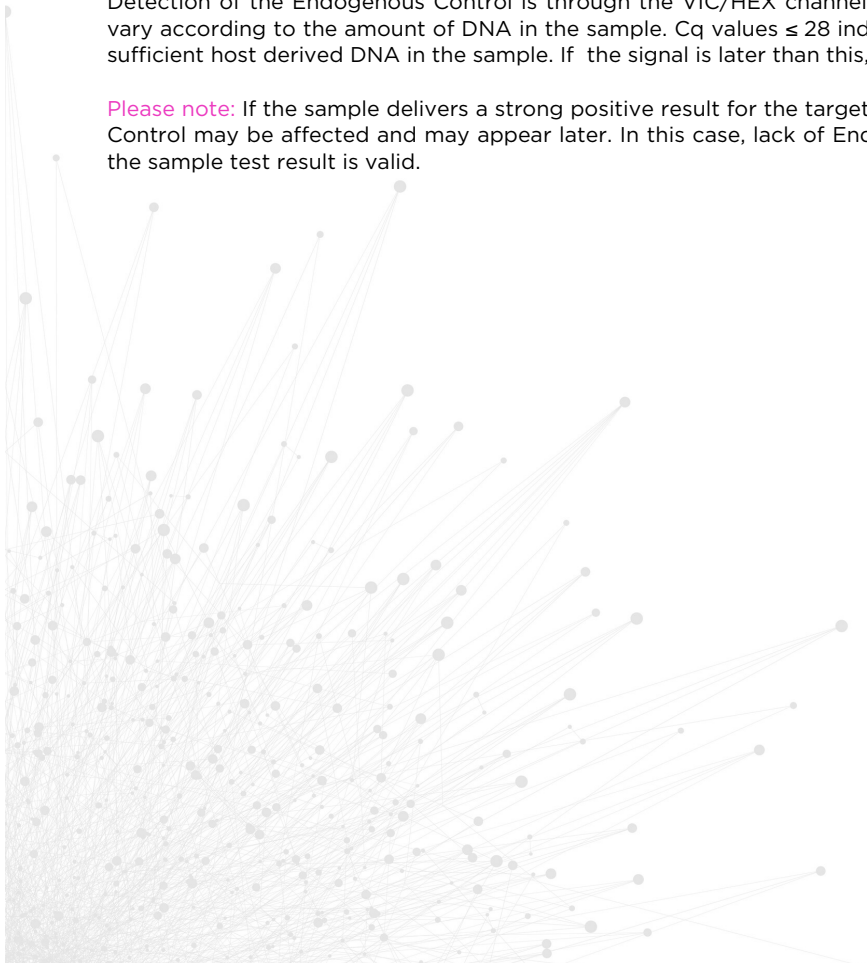



## PRODUCT SPECIFICATIONS

### Storage

Store at -20°C from arrival. The product's shelf life is outlined as an expiry date on the pouch label.

### Suitable input material

This qPCR test will work well with any source of good quality nucleic acid. Good quality is defined as nucleic acid with high integrity (not degraded). Poor quality input nucleic acid is a leading cause of test failure.

### Regulatory status

This product has been developed for Research Use Only (RUO) and is not intended for diagnostic use. It should not be used for diagnosis of disease or infection unless specifically approved by the regulatory authorities in the country of use.

### Quality Control

In accordance with the YouSeq Ltd ISO EN 13485-certified Quality Management System, each lot of *Corynebacterium minutissimum* qPCR Test Kit is tested against predetermined specifications to ensure consistent product quality. The primers/probe(s) typically demonstrate  $\geq 95\%$  in silico specificity to their intended target and periodically checked against newly available sequence information to maintain their detection profile.

### Technical Assistance

For customer support, please contact:

e-mail: [support@youseq.com](mailto:support@youseq.com)

phone: +44 (0)333 577 6697

### Trademarks and Disclaimers

YouSeq®, Tetra™, ROX™, VIC™.

Registered names, trademarks, etc. used in this document, even if not specifically marked as such, are not to be considered unprotected by law.

Not available in all countries.

© 2025 YouSeq Ltd; all rights reserved.

